

PD-13
Pediatric Plates SS
2.7 mm Set

Hip Plates SS 2.7 mm

The Pediatric Hip Plates are designed for stable fixation of varus, valgus or rotational osteotomies and trauma applications in pediatric orthopedics and to meet the specific requirements of pediatric orthopedic surgery. The Pediatric Hip Plates have a universal design for the left and right femur. The head of the plate features threaded holes for locking screws that either angle into the femoral neck or angle parallels to the growth plate in the distal femur in place of the traditional blade. In the proximal femur plates, an additional, diverging calcar screw ensures increased fixation in the bone. The 100°, 110°, and 120° plates are designed with an offset for osteotomies.

Indications: This plate is intended for use in infants up to 3 years, depending on body weight and bone quality. Neglected dislocation of the hip in combination with open reduction. Developmental coxa valga. Severe hip dysplasia.

Plate Profile: 2.7 mm Plate Profile: Thickness; 3.0 mm, Width: 8.0 mm, Distance of holes: 7.5 mm. The hexalobular recess screws are used with the plates.

TST Code	Barcode	Description	Angle	Offset	Material	Length(mm)
31610310027	8680858441204	LC PEDIATRIC HIP PLATE S.S. 2,7 MM	100° SHORT	6	SS	50
31610311027	8680858441211	LC PEDIATRIC HIP PLATE S.S. 2,7 MM	110° SHORT	6	SS	50
31610312427	8680858441235	LC PEDIATRIC HIP PLATE S.S. 2,7 MM	120° TRAUMA 4H	-	SS	55



Hip Plates SS 2.7 mm



Trauma Plate 120°



Varus 100°



Varus 110°



Trauma Plate 120°

Implants & Instruments

No	TST Code	Barcode	Description	Qty
1	00000307100	8680858444854	PEDIATRIC 2.7 MM PLATE 1.DESIGN TRAY AL	1
2	31610310027	8680858441204	LC PEDIATRIC HIP PLATE STEEL 2.7 MM X 100°	2
3	31610311027	8680858441211	LC PEDIATRIC HIP PLATE STEEL 2.7 MM X 110°	2
4	31610312427	8680858441235	LC PEDIATRIC HIP PLATE STEEL 2.7 MM X 120° STD 4 H	2
5	00000202700	8680858425761	PEDIATRIC PLATE 2.7 MM SCREW BOX	1
6	02005000500	8698673497023	DEPTH GAUGE (0-50 MM)	1
7	02025010323	8699931010947	SCREWDRIVER QUICK TIP T8X80MM	2
8	02005100020	8699931019094	DRILL GUIDE 2 MM (HANDLED)	1
9	02005800020	8699931029581	THREADED DRILL SLEEVE Ø 2X40 MM FOR SMALL HEAD SCR	2
10	02010101001	8698673441187	SOFT SCREW DRIVER QUICK SMALL	1
11	03160200015	8699931003901	PLATE GUIDE BLOCK HOLDER	1
12	02055027000	8680858425747	LENGTH MEASUR. DEV. - LC PEDIATRIC HIP PLTS. 2.7	1
13	04551000160	8699931028805	K-WIRE TUBE Ø10XØ8X160 MM	1
	22310130020	8698673450400	BONE DRILL 2.0X130 MM	2
	23412150020	8698673452701	KIRSCHNER WIRE TREADED POINT 2X150 MM	2
	23410150120	8698673453104	KIRSCHNER WIRE TROCAR POINT 2X150 MM	2



LOCKING SCREW T8 SS 2.7 MM

TST Code	Barcode	Length	Qty
21617120827	8680858445646	12	8
21617140827	8680858445653	14	8
21617160827	8680858445660	16	8
21617180827	8680858445677	18	8
21617200827	8680858445684	20	8
21617220827	8680858445691	22	8
21617240827	8680858445707	24	8
21617260827	8680858445714	26	8
21617280827	8680858445721	28	8
21617300827	8680858445738	30	4
21617320827	8680858445745	32	4
21617340827	8680858445752	34	4
21617360827	8680858445769	36	4
21617380827	8680858445776	38	4
21617400827	8680858445783	40	4





SCREW T8 SS 2.7 MM

TST Code	Barcode	Length	Qty
21610120827	8680858445790	12	4
21610140827	8680858445806	14	4
21610160827	8680858445813	16	4
21610180827	8680858445820	18	4
21610200827	8680858445837	20	4
21610220827	8680858445844	22	4
21610240827	8680858445851	24	4
21610260827	8680858445868	26	4
21610280827	8680858445875	28	4
21610300827	8680858445882	30	2
21610320827	8680858445899	32	2
21610340827	8680858445905	34	2
21610360827	8680858445912	36	2
21610380827	8680858445929	38	2
21610400827	8680858445936	40	2

No	TST Code	Barcode	Description	Qty
14	00000307200	8680858444861	PEDIATRIC 2.7 MM PLATE 2.DESIGN TRAY AL	1
15	03030000195	8698673493223	BONE-HOLDING FORCEPS 195 MM	1
16	03160000195	8699931003833	REDUCTION FORCEPS 195 MM (SERRATED)	1
17	03161010001	8699931008968	OSTEOTOMY RULER	1
18	02060000020	8699931003857	2MM K-WIRE GUIDE (ADJUSTABLE)	1
19	03160100015	8699931003871	POSITION DEVICE(GUIDE BLOCK)	1
20	03160100270	8680858431960	100° & 110° & 120° & 150° PLATE 2.7 MM GUIDE BLOCK	1
21	03381016030	8699931023800	CHISEL BLADE 16X30 MM STRAIGHT	1
	03381011030	8699931023794	CHISEL BLADE 11X30 MM STRAIGHT	1
22	03381000001	8699931023770	CHISEL HANDLE	1
23	04015000025	8699931010725	ALLEN WRENCH Ø 2.5 MM	1
	00270270110	8699931010763	CONTAINER 270X270X110 MM	1

Implants & Instruments

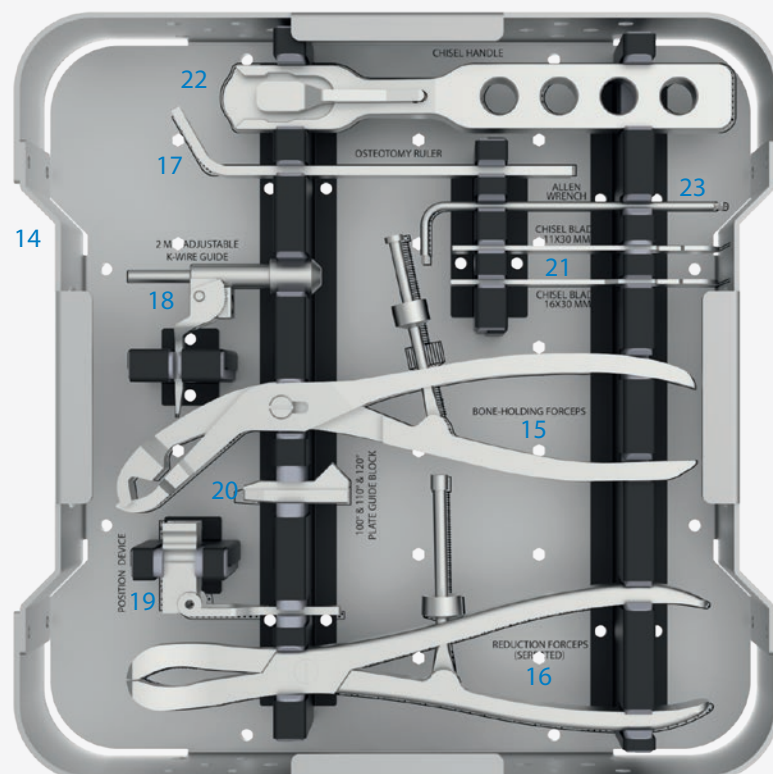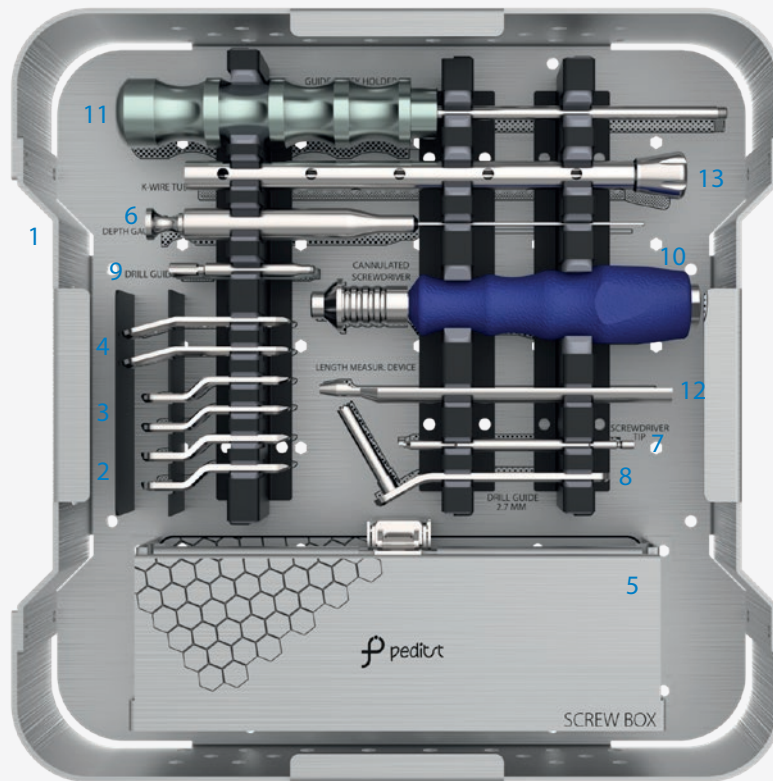

PD-08
Pediatric Plates Ti
2.7 mm Set

HIP PLATES 2.7 mm

The Pediatric Hip Plates are designed for stable fixation of varus, valgus or rotational osteotomies and trauma applications in pediatric orthopedics and to meet the specific requirements of pediatric orthopedic surgery. The Pediatric Hip Plates have a universal design for the left and right femur. The head of the plate features threaded holes for locking screws that either angle into the femoral neck or angle parallels to the growth plate in the distal femur in place of the traditional blade. In the proximal femur plates, an additional, diverging calcar screw ensures increased fixation in the bone. The 100°, 110°, and 120° plates are designed with an offset for osteotomies.


Indications: This is intended for use in infants up to 3 years, depending on body weight and bone quality. Neglected dislocation of the hip in combination with open reduction. Developmental coxa valga. Severe hip dysplasia.

Plate Profile: Thickness: 3.0 mm, Width: 8.0 mm, Distance of holes: 7.5 mm.


TST Code	Barcode	Description	Angle	Material	Offset	Length(mm)
31620310027	8699931027136	LC PEDIATRIC HIP PLATE TI 2.7 MM	100° SHORT	Ti	6	50
31620311027	8699931027143	LC PEDIATRIC HIP PLATE TI 2.7 MM	110° SHORT	Ti	6	50
31620312427	8699931027167	LC PEDIATRIC HIP PLATE TI 2.7 MM	120° TRAUMA	Ti	-	55



HIP PLATES 2.7 mm

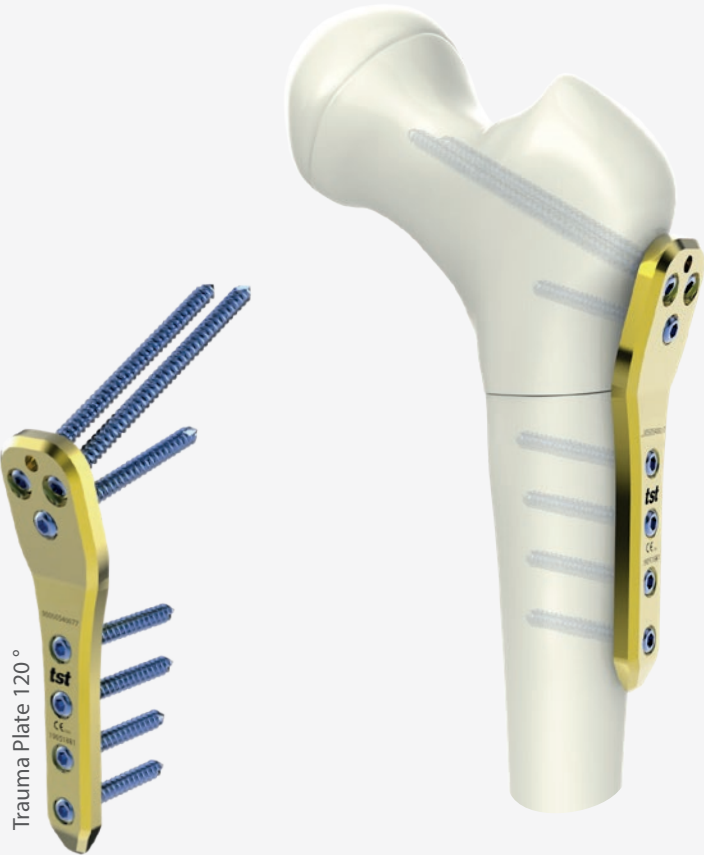

Trauma Plate 120°

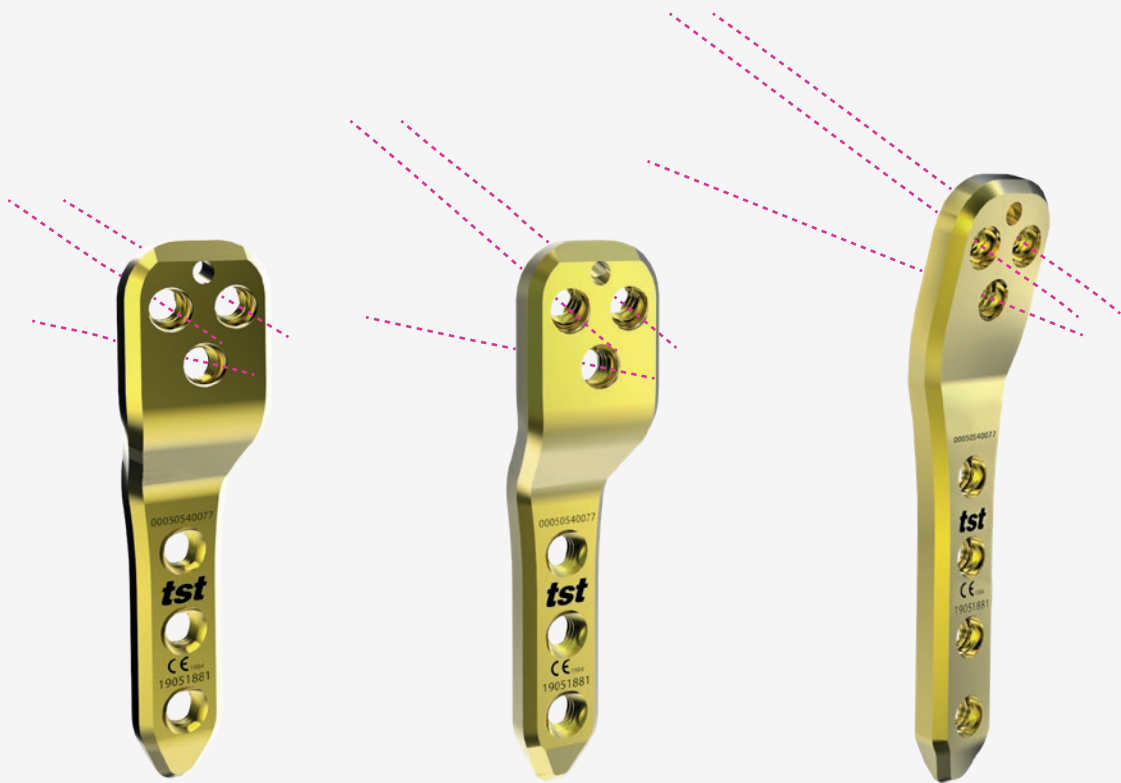

Varus 100°

Varus 110°

Trauma Plate 120°

Implants & Instruments

No	TST Code	Barcode	Description	Qty
1	00000307100	8680858444854	PEDIATRIC 2.7 MM PLATE 1.DESIGN TRAY AL	1
2	31620310027	8699931027136	LC PEDIATRIC HIP PLATE TI 2.7 MM X 100° SHORT	2
3	31620311027	8699931027143	LC PEDIATRIC HIP PLATE TI 2.7 MM X 110° SHORT	2
4	31620312427	8699931027167	LC PEDIATRIC HIP PLATE TI 2.7 MM X 120° TRAUMA	2
5	00000202700	8680858425761	PEDIATRIC PLATE 2.7 MM SCREW BOX	1
6	02005000500	8698673497023	DEPTH GAUGE (0-50 MM)	1
7	02000480017	8699931029253	SCREW DRIVER QUICK HEX. TIP Ø 1.7X80 MM	2
8	02005100020	8699931019094	DRILL GUIDE 2 MM (HANDLED)	1
9	02005800020	8699931029581	THREADED DRILL SLEEVE Ø 2X40 MM FOR SMALL HEAD SCR	2
10	02010101001	8698673441187	SOFT SCREW DRIVER QUICK SMALL	1
11	03160200015	8699931003901	PLATE GUIDE BLOCK HOLDER	1
12	02055027000	8680858425747	LENGTH MEASUR. DEV. - LC PEDIATRIC HIP PLTS. 2.7	1
13	04551000160	8699931028805	K-WIRE TUBE Ø10XØ8X160 MM	1
	22310130020	8698673450400	BONE DRILL 2.0X130 MM	2
	23412150020	8698673452701	KIRSCHNER WIRE TREADED POINT 2X150 MM	2
	23410150120	8698673453104	KIRSCHNER WIRE TROCAR POINT 2X150 MM	2



SMALL HEAD LOCK. SCREW TI 2.7 MM

TST Code	Barcode	Length	Qty
24627100027	8699931015539	8	8
24627120027	8699931007220	8	8
24627140027	8699931007213	8	8
24627160027	8699931007206	8	8
24627180027	8699931007190	8	8
24627200027	8699931007183	8	8
24627220027	8699931007176	8	8
24627240027	8699931005912	8	8
24627260027	8699931005929	8	8
24627280027	8699931005936	8	8
24627300027	8699931005943	4	4
24627320027	8699931005950	4	4
24627340027	8699931005967	4	4
24627360027	8699931005974	4	4
24627380027	8699931005981	4	4
24627400027	8699931005998	4	4



SMALL HEAD SCREW TI 2.7 MM

TST Code	Barcode	Length	Qty
24620100027	8699931015140	4	4
24620120027	8699931006841	4	4
24620140027	8699931006834	4	4
24620160027	8699931006827	4	4
24620180027	8699931006810	4	4
24620200027	8699931006803	4	4
24620220027	8699931006797	4	4
24620240027	8699931006162	4	4
24620260027	8699931006155	4	4
24620280027	8699931006148	4	4
24620300027	8699931006131	2	2
24620320027	8699931006124	2	2
24620340027	8699931006117	2	2
24620360027	8699931006100	2	2
24620380027	8699931006094	2	2
24620400027	8699931006087	2	2

No	TST Code	Barcode	Description	Qty
14	00000307200	8680858444861	PEDIATRIC 2.7 MM PLATE 2.DESIGN TRAY AL	1
15	03030000195	8698673493223	BONE-HOLDING FORCEPS 195 MM	1
16	03160000195	8699931003833	REDUCTION FORCEPS 195 MM (SERRATED)	1
17	03161010001	8699931008968	OSTEOTOMY RULER	1
18	02060000020	8699931003857	2MM K-WIRE GUIDE (ADJUSTABLE)	1
19	03160100015	8699931003871	POSITION DEVICE(GUIDE BLOCK)	1
20	03160100270	8680858431960	100° & 110° & 120° & 150° PLATE 2.7 MM GUIDE BLOCK	1
21	03381016030	8699931023800	CHISEL BLADE 16X30 MM STRAIGHT	1
	03381011030	8699931023794	CHISEL BLADE 11X30 MM STRAIGHT	1
22	03381000001	8699931023770	CHISEL HANDLE	1
23	04015000025	8699931010725	ALLEN WRENCH Ø 2.5 MM	1
	00270270110	8699931010763	CONTAINER 270X270X110 MM	1

Implants & Instruments

