

PD-10
BUTTERFLY GROWTH
PLATES SET




BUTTERFLY
GROWTH

BUTTERFLY GROWTH PLATES


Butterfly Growth Plate is available in two sizes, 13 and 16 mm depending on the distance between the screw holes. The system is composed of biocompatible Titanium alloy. 16 mm plate with 2- and 4-holes options are available. Locking cannulated screws and locking peg screws are both included in the set. Butterfly Growth Plate has a durable structure against to the force for hemi-epiphyseodesis applications. It can be applied with a precise and simple surgical technique which provides an easy adaptation to different bone morphologies. Each hinge of plates allows up to 18 degrees of angulation. The hinged structure of the system ensures that the plate is fully aligned with the bone surface, providing an ideal positioning on the growth plate. Butterfly Growth Plate facilitates applications by providing the desired full anatomical fit to inclined surfaces. The pull-out risk of the PEG screw is the same as for the threaded screw. They also have minimal damage to the bone stock.

Indications: Butterfly Growth Plate has many pathological indications such as posttraumatic, metabolic, neuromuscular, idiopathic, and congenital pathologies. It is preferred for long bone deformities and corrections of unequal legs.

TST Code	Barcode	Description	Material	Size
33620000122	8680858426768	BUTTERFLY GROWTH PLATE 2 HOLES 13 MM	Ti	13 mm
33620000162	8680858426775	BUTTERFLY GROWTH PLATE 2 HOLES 16 MM	Ti	16 mm
33620000164	8680858426805	BUTTERFLY GROWTH PLATE 4 HOLES 16 MM	Ti	16 mm

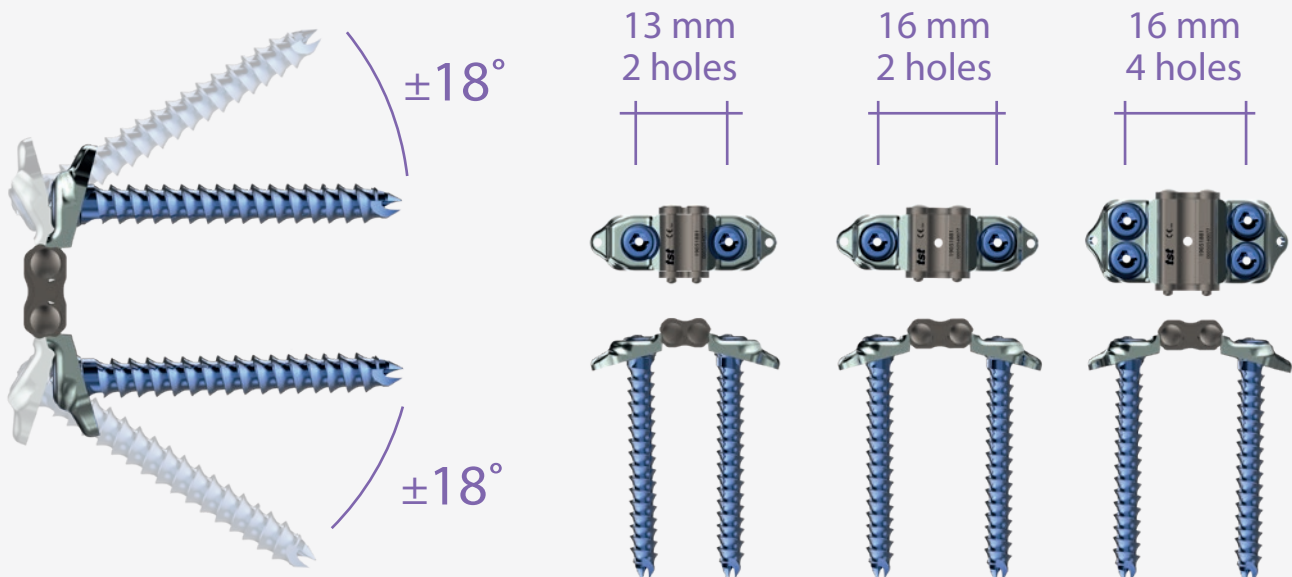

BUTTERFLY GROWTH PLATES

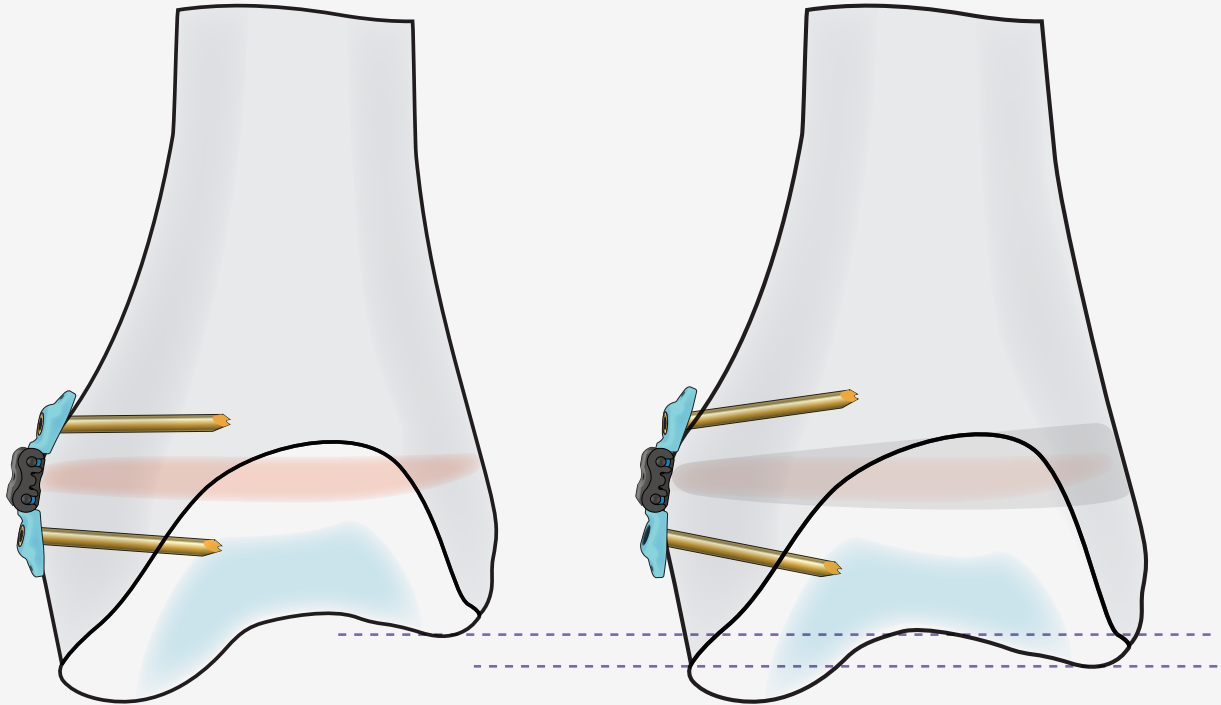

LOCK. CANN. PEG SCREW Ti Ø 3.5 MM

TST Code	Barcode	Length	Qty
20127002035	8680858427772	20 mm	2
20127002235	8680858427789	22 mm	2
20127002435	8680858427796	24 mm	2
20127002635	8680858427802	26 mm	2
20127002835	8680858427819	28 mm	2
20127003035	8680858427826	30 mm	2
20127003535	8680858427833	35 mm	2
20127004035	8680858427840	40 mm	2



LOCK. CANNULATED SCREW Ti Ø 4 MM

TST Code	Barcode	Length	Qty
20427002040	8680858413737	20 mm	2
20427002240	8680858413744	22 mm	2
20427002440	8680858413751	24 mm	2
20427002640	8680858413768	26 mm	2
20427002840	8680858413775	28 mm	2
20427003040	8680858413782	30 mm	2
20427003540	8680858413799	35 mm	2
20427004040	8680858413805	40 mm	2

Implants & Instruments

No	TST Code	Barcode	Description	Qty
1	00000304100	8680858437771	DESIGN TRAY 250X250 MM	1
2	00000304010	8680858437764	BUTTERFLY GROWTH PLATE&SCREW BOX	1
3	33620000122	8680858426768	BUTTERFLY GROWTH PLATE 2 HOLES (TI) 13 MM	4
4	33620000162	8680858426775	BUTTERFLY GROWTH PLATE 2 HOLES (TI) 16 MM	4
5	33620000164	8680858426805	BUTTERFLY GROWTH PLATE 4 HOLES (TI) 16 MM	4
6	03362020100	8680858437757	PENSET FOR BUTTERFLY GROWTH PLATE	1
7	3362010010	8680858437733	DRILL SLEEVE HOLDER	1
8	03362010100	8680858437740	SCREW EXTRACTOR FOR BUTTERFLY GROWTH PLT	1
9	03362006532	8680858437719	THREADED DRILL SLEEVE 3.2X65 MM	2
10	03362007512	8680858437726	K-WIRE SLEEVE 1.2X100 MM	2
11	01194003009	8699931028157	BONE HAMMER SMALL WITH SLOTTED END	1
12	3356200100	8680858442256	LENGTH MEASURING DEVICE FOR K-WIRE 1.2X150 MM	1
13	02010101001	8698673441187	SOFT SCREW DRIVER QUICK SMALL	1
14	02050100425	8699931005417	CANNULATED SCREW DRIVER QUICK TIP 2.5XØ5X130MM	2
15	04551208190	8699931029550	K-WIRE TUBE Ø12XØ8X190 MM	1
	23410150112	8698673458840	KIRSCHNER WIRE TROCAR POINT 1.2X150 MM	4
	23212530170	8680858426355	CANNULATED DRILL 3.2X1.25X30X170 MM	2
	00270270080	8699931009019	CONTAINER 270X270X80 MM	1

Implants & Instruments



PD-10 Butterfly Growth Plates Set

