

EF-18
PEDIATRIC DAF
FIXATORS



DAF FIXATOR

It is a single plan external fixation tool. The ball-joint structure connects the clamp to the fixator body as an articular structure and allows the final reduction after the implantation of the fixator. The Fixators Schanz screw holder clamps have five slots for the Schanz screws. It is primarily used in acute trauma cases. In some cases, if it is needed, it allows joint mobilization. The system allows correction of postoperative reduction. Fixators versatility allows simple angular adjustments. It provides an opportunity to move 30° in all planes. A full turn of the DC unit allows 1 mm movement.

Indications:

Treatment of the fractures of femur and tibia and deformity corrections in pediatric patients.

TST Code	Barcode	Description	Type	Length (mm)	Lengthning (mm)
40840010101	8698673456525	DYN. MULTIAXIAL FIX.	Standard	175	25
40840020101	8698673456518	DYN. MULTIAXIAL FIX.	Long	200	47



DAF WITH 135° CLAMP FIXATOR

DAF WITH 135° CLAMP FIXATOR has Schanz screw holes opened in a 135 ° angle to the long axis of the clamp. It allows the placement of the Schanz screws to the Femur head when the fixator body is positioned parallel to the long axis of Femur. The 135 ° clamp has four screw holes.



TST Code	Barcode	Description
40840000300	8698673456709	135° CLAMP

TST Code	Barcode	Description
40840010151	8698673457980	DYN. MULTIAXIAL WITH 135° CLAMP FIX.
40840020151	8698673457973	DYN. MULTIAXIAL WITH 135° CLAMP FIX.

Type	Length (mm)	Lengthning (mm)
Standard	175	25
Long	200	47

DAF WITH T CLAMP FIXATOR

T-clamp with the ball-joint structure could be used instead of the straight clamp at the dynamic multi-axial fixator. T-Clamp gives the opportunity to place the Schanz screws in a uni-planar configuration to epiphyseal and metaphyseal region. It has a design that 5 screws could be applied horizontally to the metaphyseal region of the bone. It is used for distal/proximal femur and distal/proximal tibia intra-articular and periarticular fractures which have enough space for Schanz screw applications.



TST Code	Barcode	Description
40840000200	8698673456716	T CLAMP

TST Code	Barcode	Description
40840010141	8698673458024	DYN. MULTIAXIAL WITH T CLAMP FIX.
40840020141	8698673458031	DYN. MULTIAXIAL WITH T CLAMP FIX.

Type	Length (mm)	Lengthning (mm)
Standard	155	25
Long	180	47

DYNAMIC ANGLED CLAMP FIXATOR

Lightweight and low-profile structured fixator. Owing to the double ball-joint structure, it provides an easy reduction. The ball joint structure can be tightened and / or loosened with a 4.0 mm L-Allen key. Suitable for extra-articular or trans articular external fixation technique. The ball-joint structure can be loosened at the trans articular application if there is a need of a Hand-wrist mobilization. It is used with 3-4 mm Schanz screws.

Indications:

It is used in the treatment of joint stiffness, articular and periarticular fractures of distal radius and deformity corrections.



TST Code	Barcode	Description	Type	Length (mm)
40640030103	8698673456570	DYNAMIC ANGLED CLAMP WRIST FIXATOR	SHORT	170
40640020102	8698673456587	DYNAMIC ANGLED CLAMP WRIST FIXATOR	STANDARD	195
40640010101	8698673456563	DYNAMIC ANGLED CLAMP WRIST FIXATOR	LONG	230

The Pediatric LRS system aims to treat successfully congenital and post-traumatic deformities for children and small statured adults. Small sizes specifically developed for patient comfort.

Indications: Limb reconstruction and lengthening due to fresh fracture, nonunion with major soft tissue defect, and bone loss with shortening.



TST Code	Barcode	Description	Type	Length (mm)
44140010250	8680858407538	PEDIATRIC LRS FIX.	LONG	250
44140020200	8680858407521	PEDIATRIC LRS FIX.	STANDARD	200
44140030150	8680858407514	PEDIATRIC LRS FIX.	SHORT	150
44140000001	8680858407545	PEDIATRIC RAIL CLAMP		

Implants & Instruments

No	TST Code	Barcode	Description	Qty
1	00000703100	8699931017267	HUMERUS - WRIST DESIGN 1. TRAY	1
2	00000601300	8699931017182	3MM\4MM SCHANZ BOX	1
	20910101503	8698673455863	DAF SCHANZ PINCORTICAL THREADED SS 3X15X100	4
	20910002003	8698673455843	DAF SCHANZ PINCORTICAL THREADED SS 3X20X100	4
	20910002503	8698673455917	DAF SCHANZ PINCORTICAL THREADED SS 3X25X100	4
	20910601503	8698673455832	DAF SCHANZ PINCORTICAL THREADED SS 3X15X60	4
	20910602003	8698673455887	DAF SCHANZ PINCORTICAL THREADED SS 3X20X60	4
	20910801503	8698673455849	DAF SCHANZ PINCORTICAL THREADED SS 3X15X80	4
	20910802003	8698673455894	DAF SCHANZ PINCORTICAL THREADED SS 3X20X80	4
	20910121503	8698673455870	DAF SCHANZ PINCORTICAL THREADED SS 3X15X120	4
	20910202003	8698673455900	DAF SCHANZ PINCORTICAL THREADED SS 3X20X120	4
	20910702004	8698673458604	DAF SCHANZ PINCORTICAL THREADED SS 4X20X70	4
	20910801504	8698673458116	DAF SCHANZ PINCORTICAL THREADED SS 4X15X80	4
	20910802004	8698673455962	DAF SCHANZ PINCORTICAL THREADED SS 4X20X80	4
	20910002004	8698673473317	DAF SCHANZ PINCORTICAL THREADED SS 4X20X100	4
	20910102504	8698673456013	DAF SCHANZ PINCORTICAL THREADED SS 4X25X100	4
	20910202004	8698673455986	DAF SCHANZ PINCORTICAL THREADED SS 4X20X120	4
	20910202504	8698673456020	DAF SCHANZ PINCORTICAL THREADED SS 4X25X120	4
	20910402004	8698673458383	DAF SCHANZ PINCORTICAL THREADED SS 4X20X140	4
	20910142504	8698673456037	DAF SCHANZ PINCORTICAL THREADED SS 4X25X140	4
3	02005100012	8698673490703	DRILL GUIDE (SINGLE)	1
4	02005100015	8699931017854	DRILL GUIDE DOUBLE WRIST FIXATOR	1
5	04015000300	8698673492493	3 MM UNIVERSAL SCHANZ T HANDLE	1
6	04015000400	8698673492509	4 MM UNIVERSAL SCHANZ T HANDLE	1
7	00000703200	8699931017274	HUMERUS - WRIST DESIGN 2. TRAY	1
8	40840010101	8698673456525	DYN. MULTIAXIAL HUMERUS FIXATOR STD.	1
9	40840020101	8698673456518	DYN. MULTIAXIAL HUMERUS FIXATOR LONG	1
10	40840000200	8698673456716	HUMERUS T CLAMP	*
11	40840000300	8698673456709	HUMERUS 135 CLAMP	1
12	40640010101	8698673456563	DYNAMIC ANGLED CLAMP WRIST FIXATOR LONG	2
13	40640020102	8698673456587	DYNAMIC ANGLED CLAMP WRIST FIXATOR STANDARD	2


* Titanium Schanz Pins (optional)

TST Code	Barcode	Description
20920101503	8698673457140	DAF SCHANZ PIN CORT.THRD.(TI) 3X15X100
20920102503	8698673455320	DAF SCHANZ PIN CORT.THRD.(TI) 3X25X100
20920062003	8698673455306	DAF SCHANZ PIN CORT.THRD.(TI) 3X20X 60
20920081503	8698673457133	DAF SCHANZ PIN CORT.THRD.(TI) 3X15X80
20920082003	8698673455313	DAF SCHANZ PIN CORT.THRD.(TI) 3X20X 80
20920083504	8698673458758	DAF SCHANZ PIN CORT.THRD.(TI) 4X35X80
20920102504	8698673455351	DAF SCHANZ PIN CORT.THRD.(TI) 4X25X100
20920122504	8698673455368	DAF SCHANZ PIN CORT.THRD.(TI) 4X25X120
20921402504	8699931028607	DAF SCHANZ PIN CORT.THRD.(TI) 4X25X140

* optional

please contact with our sales specialist for more information

Implants & Instruments

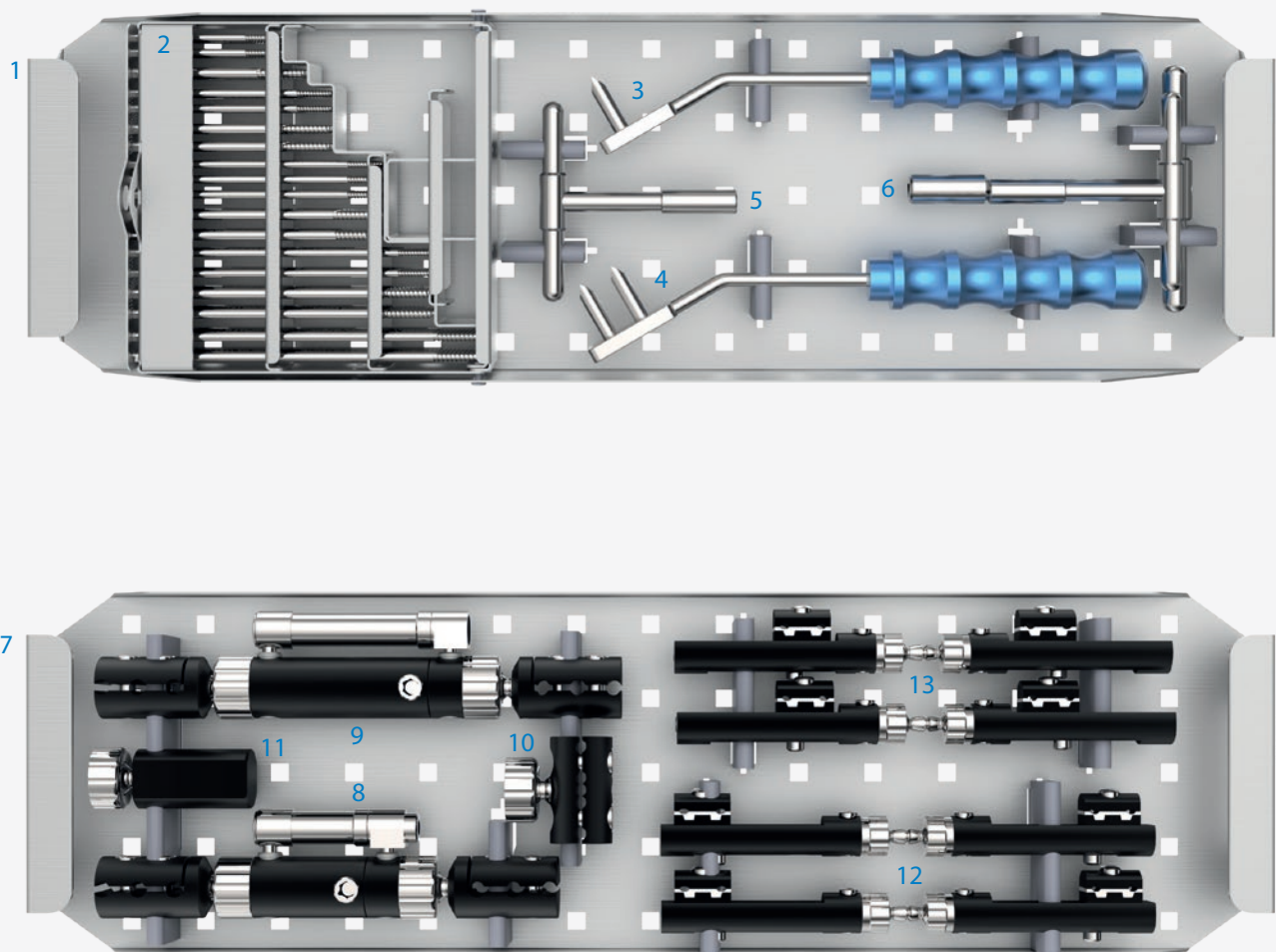

EF-18 DAF & Wrist Fixator Set



Implants & Instruments

No	TST Code	Barcode	Description	Qty
14	00000703300	8699931017281	HUMERUS - WRIST DESIGN 3. TRAY	1
15	02005100025	8698673490178	AWL BLUNT TIPPED	1
16	02005100023	8698673490185	AWL TROCAR TIPPED	1
	04100020044	8680858431984	PROTECTIVE CAP FOR 4 MM SCHANZ	10
	04100020023	8680858431977	PROTECTIVE CAP FOR 2-3 MM SCHANZ	10
17	02006000200	8699931015454	DRILL GUIDE SLEEVE 2 MM	1
18	02006000250	8699931017885	DRILL GUIDE SLEEVE 2.5MM	1
19	02006000320	8699931023411	DRILL GUIDE SLEEVE 3.2 MM	1
20	02005120000	8698673490833	DRILL BIT WITH STOP Ø 2 MM	2
21	02005125000	8698673492950	DRILL BIT WITH STOP Ø 2.5 MM	2
22	02005132000	8698673495432	DRILL BIT WITH STOP Ø 3.2 MM	2
23	04015000020	8699931015430	ALLEN WRENCH Ø 2 MM	1
24	04015000040	8698673490055	ALLEN WRENCH Ø 4 MM	2
25	04015000050	8698673496163	ALLEN WRENCH Ø 5 MM	1
	00500160150	8699931017199	CONTAINER 500X160X150 MM	1
* 26	00000703400	8680858451777	PEDI. LRS TRAY	1
27	44140010250	8680858407538	PEDIATRIC LRS FIX. LONG 250 MM	1
28	44140020200	8680858407521	PEDIATRIC LRS FIX. STANDARD 200 MM	1
29	44140030150	8680858407514	PEDIATRIC LRS FIX. SHORT 150 MM	1
30	44140000001	8680858407545	PEDIATRIC RAIL CLAMP	1

* optional

please contact with our sales specialist for more information

Implants & Instruments

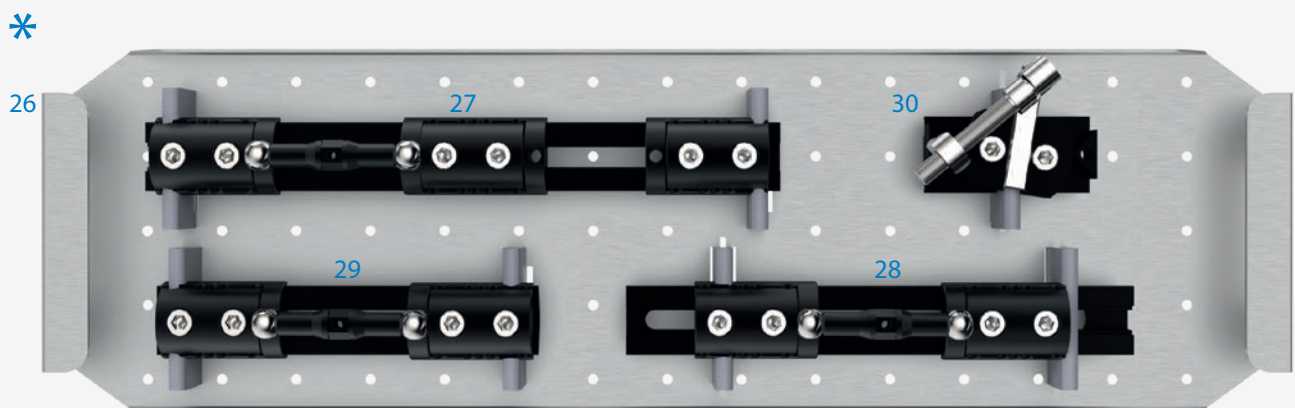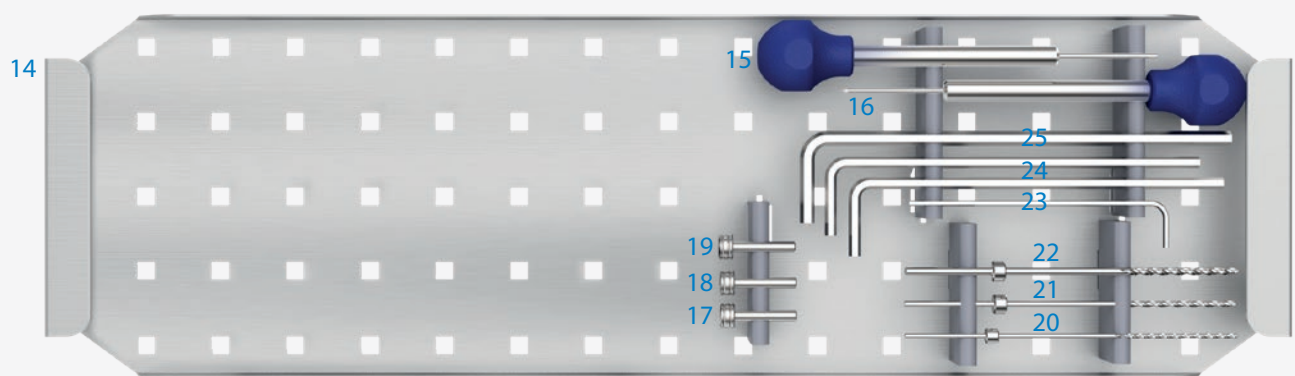

EF-18 DAF & Wrist Fixator Set

