

PD-09
DERONA 360°
NAIL SET

DeRoNa 360°

DERONA 360° NAIL

Indications: DeRoNa 360 is designed for proximal femur derotation osteotomy treatments, proximal region fractures, malunion, nonunion and deformities in the child age.

Features:

*DeRoNa 360° is applied from the lateral trochanteric entry point.

*DeRoNa 360° has 12° of valgus angle. The length of the proximal slope is 32 mm.

*DeRoNa 360° consists of proximal and distal parts. Proximal slope of the nail is designed to make the easiest entry in the sagittal and coronal plane.

*Proximal and distal components are loosened for controlled rotational correction to provide the proximal femur anteversion correction. Distal and Proximal components are connected with the connection bolt to form the rigid nail structure.

*The square-shaped proximal part is designed to increase the rotational stability in the proximal femur.

*The proximal part of DeRoNa 360° has three different width dimensions (10-11-12 mm).

*The locking screws can be applied from 3 different positions of the proximal part.


*The locking screws are available in 20 different sizes from 20 to 70 mm.

*The distal component is available in three sizes (100-125-150 mm) and three different diameters (Ø 8-9-10mm).

*DeRoNa 360° has 1 piece on a 100 mm long distal component, 2 pieces on 125 mm, and 3 static locking holes on the 150 mm long component.

*End Caps are available in four different sizes as 0-5-10-15 mm.

*The surgical instruments (nail guides, reamers, inserters, etc.) are placed in their trays in a container suitable for sterilization.

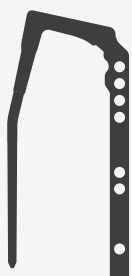

guide





rotation



locking system



DERONA 360° NAIL

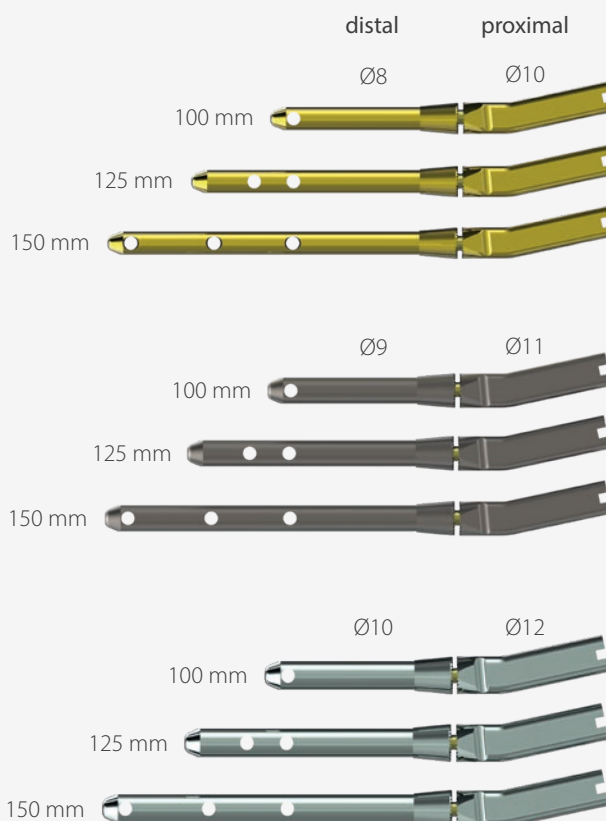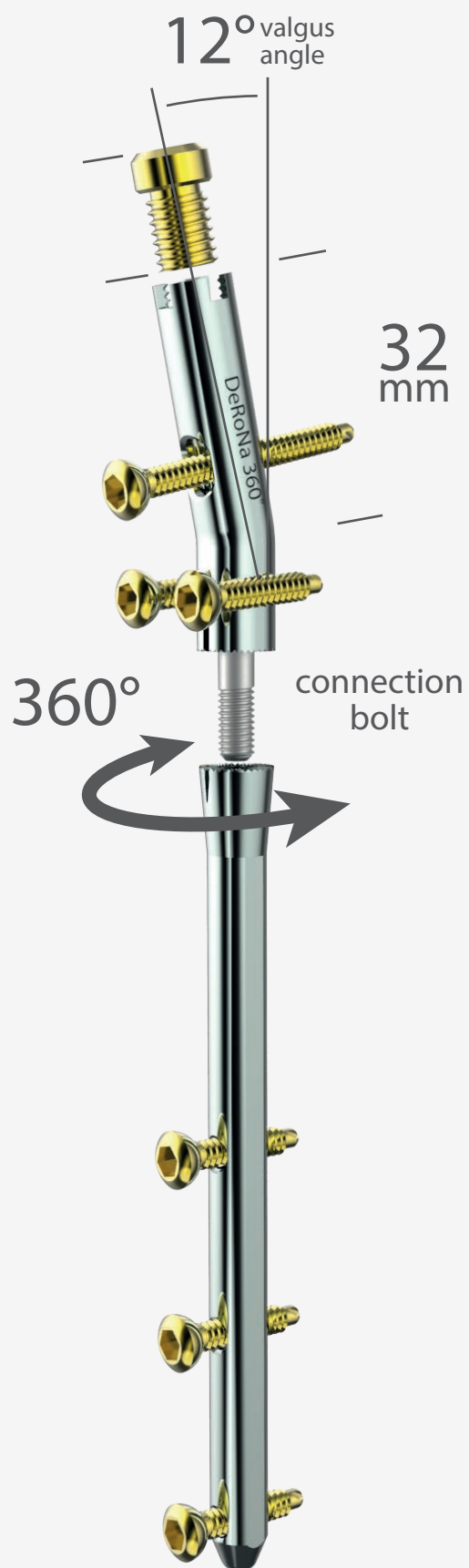



Cortex Screw for Nails Ti Ø 4.0 mm



Implants & Instruments

No	TST Code	Barcode	Description	Qty
1	00000802210	8680858437160	DERONA IMPLANT TRAY	1
2	80220210010	8680858404575	DERONA 360 8x100 MM	2
	80220220010	8680858404582	DERONA 360 8x125 MM	2
	80220250010	8680858404599	DERONA 360 8x150 MM	2
3	80220210011	8680858404605	DERONA 360 9x100 MM	2
	80220220011	8680858404612	DERONA 360 9x125 MM	2
	80220250011	8680858404629	DERONA 360 9x150 MM	2
4	80220220012	8680858404636	DERONA 360 10x100 MM	2
	80220250012	8680858404643	DERONA 360 10x125 MM	2
	80220000000	8680858408122	DERONA 360 10x150 MM	2
5	83021800053	8699931026900	DERONA NECK LOCK IMPLANT 11 MM	2
6	83021000053	8699931026917	DERONA NECK LOCK IMPLANT 12 MM	2
7	83021600053	8699931026894	DERONA NECK LOCK IMPLANT 10 MM	2
8	00000802220	8680858437153	DERONA SCREW BOX	1
	22124040018	8680858422340	CORTEX SCREW FOR NAILS TI Ø4.0X18 MM	3
	22124040020	8680858422357	CORTEX SCREW FOR NAILS TI Ø4.0X20 MM	3
	22124040022	8680858422364	CORTEX SCREW FOR NAILS TI Ø4.0X22 MM	3
	22124040024	8680858422371	CORTEX SCREW FOR NAILS TI Ø4.0X24 MM	3
	22124040026	8680858422388	CORTEX SCREW FOR NAILS TI Ø4.0X26 MM	3
	22124040028	8680858422395	CORTEX SCREW FOR NAILS TI Ø4.0X28 MM	3
	22124040030	8680858422401	CORTEX SCREW FOR NAILS TI Ø4.0X30 MM	3
	22124040032	8680858422418	CORTEX SCREW FOR NAILS TI Ø4.0X32 MM	3
	22124040034	8680858422425	CORTEX SCREW FOR NAILS TI Ø4.0X34 MM	3
	22124040036	8680858422432	CORTEX SCREW FOR NAILS TI Ø4.0X36 MM	3
	22124040038	8680858422449	CORTEX SCREW FOR NAILS TI Ø4.0X38 MM	3
	22124040040	8680858422456	CORTEX SCREW FOR NAILS TI Ø4.0X40 MM	3
	22124040044	8680858422470	CORTEX SCREW FOR NAILS TI Ø4.0X44 MM	3
	22124040048	8680858422494	CORTEX SCREW FOR NAILS TI Ø4.0X48 MM	3
	22124040052	8680858422517	CORTEX SCREW FOR NAILS TI Ø4.0X52 MM	3
	22124040056	8680858422531	CORTEX SCREW FOR NAILS TI Ø4.0X56 MM	3
	22124040060	8680858422555	CORTEX SCREW FOR NAILS TI Ø4.0X60 MM	3
	22124040064	8680858422579	CORTEX SCREW FOR NAILS TI Ø4.0X64 MM	3
	22124040068	8680858422593	CORTEX SCREW FOR NAILS TI Ø4.0X68 MM	3
	22124040070	8680858422609	CORTEX SCREW FOR NAILS TI Ø4.0X70 MM	3
9	00000209001	8699931028409	END CAP SCREW BOX	1
	81520010011	8680858408221	DERONA CONNECTION BOLT	6
	85020000000	8680858406753	DERONA END CAP 0 MM	2
	86020000005	8680858402243	DERONA END CAP 5 MM	2
	86020000010	8680858402250	DERONA END CAP 10 MM	2
	86020000015	8680858402267	DERONA END CAP 15 MM	2
10	00000802230	8680858437221	DERONA 1.INSTRUMENT TRAY	1
11	08022001065	8680858437177	REAMER 6.5X10X220 MM	1
12	08022001165	8680858437184	REAMER 6.5X11X220 MM	1
13	08022001265	8680858437191	REAMER 6.5X12X220 MM	1
14	08022011012	8680858437207	PROTECTION SLEEVE 12 MM	1
15	08022000060	8680858437214	K-WIRE GUIDE Ø2X12 MM	1
16	01610000105	8699931017113	FLEXIBLE REAMER Ø 10.5	1
17	08050000115	8699931017120	FLEXIBLE REAMER Ø 11.5	1
18	08050000125	8699931017137	FLEXIBLE REAMER Ø 12.5	1
19	01193000023	8698673493780	T QUICK HANDLE	1

Implants & Instruments

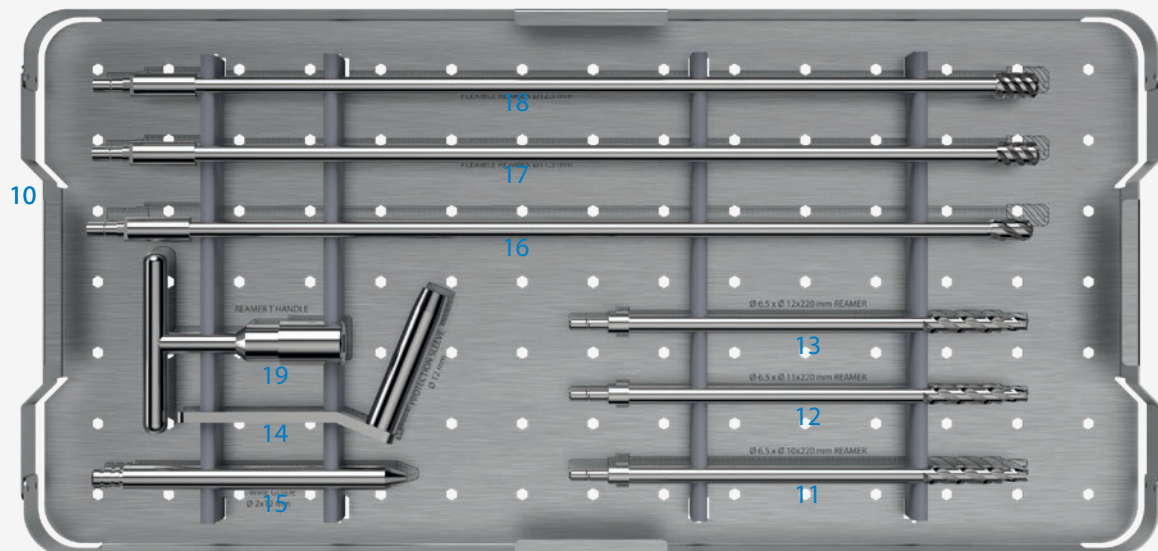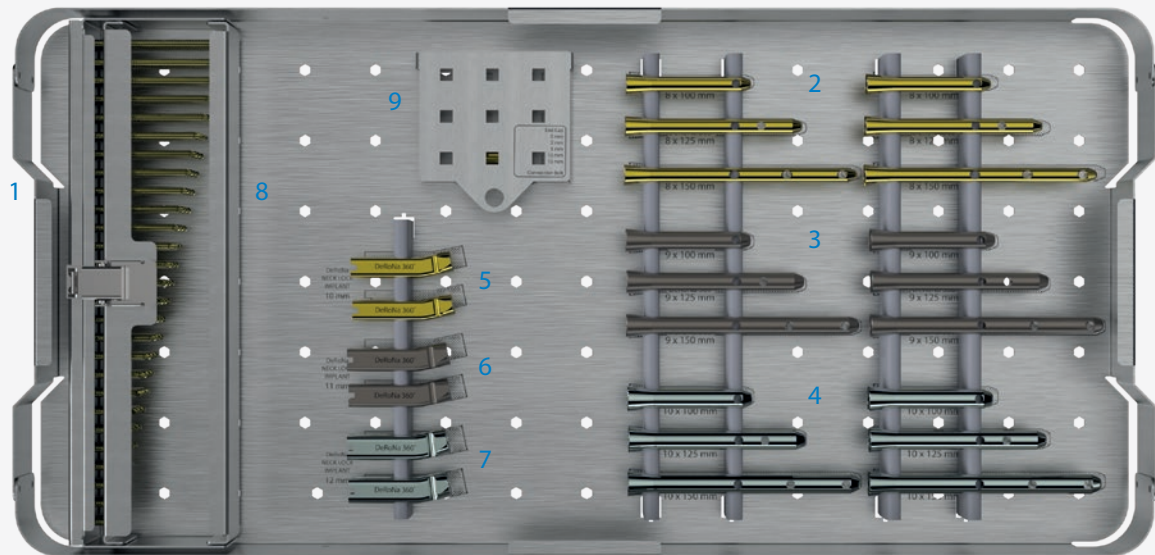

Implants & Instruments

No	TST Code	Barcode	Description	Qty
20	00000802240	8680858437320	DERONA 2.INSTRUMENT TRAY	1
21	08022001000	8680858437238	DERONA INSERTION HANDLE	1
22	08022001020	8680858437245	DERONA NAIL HOLDER CONNECTION SCREW	2
23	04015000050	8698673496163	ALLEN WRENCH Ø 5 MM	2
24	02000408150	8680858440528	SCREW DRIVER QUICK HEX. TIP 4X8X150 MM	2
25	02060406245	8680858437641	CANNULATED SCREWDRIVER Ø 4.0 MM	1
26	08022015035	8680858437252	SCREWDRIVER CROSS TIP 3.0 MM	2
27	02010101002	8698673493308	SOFT SCREW DRIVER QUICK LARGE	1
28	01191000040	8699931008739	SCREW DRIVER TORQUE 3NM	1
29	08022000020	8680858437269	DERONA SCREW SLEEVE	1
30	08022000030	8680858437276	DERONA DRILL GUIDE	1
31	08022200080	8680858437283	DERONA TROCAR 8 MM	1
32	08022000050	8680858437290	ALIGNMENT PIN	2
33	08022002508	8680858437306	DEROTATION MEASURER	2
34	08022000040	8680858437313	DEROTATION MEASURER GUIDE	1
35	00000802250	8680858437344	DERONA 3.INSTRUMENT TRAY	1
36	04551080310	8699931030815	K-WIRE TUBE Ø10XØ8X310 MM	2
	23412300020	8698673452732	KIRSCHNER WIRE TREADED POINT 2X300 MM	2
	01610030032	8680858424184	GRADUATED DRILL BIT Ø3.2X300 MM	2
37	04551208160	8699931028720	K-WIRE TUBE Ø12XØ8X160 MM	1
	23410120150	8698673452978	KIRSCHNER WIRE TROCAR POINT 1.5X120 MM	4
38	02001400100	8680858430307	DEPTH GAUGE - HIN PROX. SCREW LENGTH	1
39	03310000001	8680858400201	T-HANDLE JACOBS CHUCK-KEY	1
40	08050000200	8680858409181	T-HANDLE MALE NAIL INSERTER	1
41	01193002009	8699931028126	BONE HAMMER - MEDIUM	1
42	01610003000	8680858424054	SURGICAL RULER S.S 300 MM	1
43	08044000012	8699931015744	WRENCH 12 MM	1
44	08022001030	8680858437337	DERONA INSERTION KNOB	1
45	08061000006	8699931017335	FEMUR INTRAMEDULLARY NAIL EXTRACTOR	1
46	01195001009	8699931028195	HINGED SLOTTED HAMMER LARGE	1
	00260250600	8699931030716	PFN KIRSCHNER WIRE GUIDE Ø2.5X600	1
	00560270125	8699931010923	CONTAINER 560X270X125 MM	2

Implants & Instruments

