


PD-03  
EIGHT GROWTH  
PLATE Set



## EIGHT GROWTH PLATES

The anatomically structured titanium alloyed plate, having two different lengths (12 and 16 mm) simplifies the fixation on the physis (growth plate) by the K-Wire holes on the centre and both ends of the plate. It provides simple, short, and not complicated application opportunities by ambulatory care method. After the operation, moving easily, normal recovering process and not restricting daily activities are possible. Flexible, non-locking screw-plate hole structure provide the screws to be separated like a joint (till 24° for each screw) without restraining the physis development. Safe deformity correction is provided with the cannulated or usual, full threaded, cancellous screws having ability to resist the existing forces here. The diameters of the screws are Ø 4 mm, the lengths are from 14 mm up to 40 mm and they are made of titanium alloy.

**Indications:** The plate is applied for gradual correction of the deformities to upper and lower extremities (at the tip of long bones) during the adolescent and infant periods (from 18 months to 17 years). Ankle valgus, elbow varus, knee valgus and varus, oblique plan flexion deformity and recurvatum cases can be evaluated in this scope.

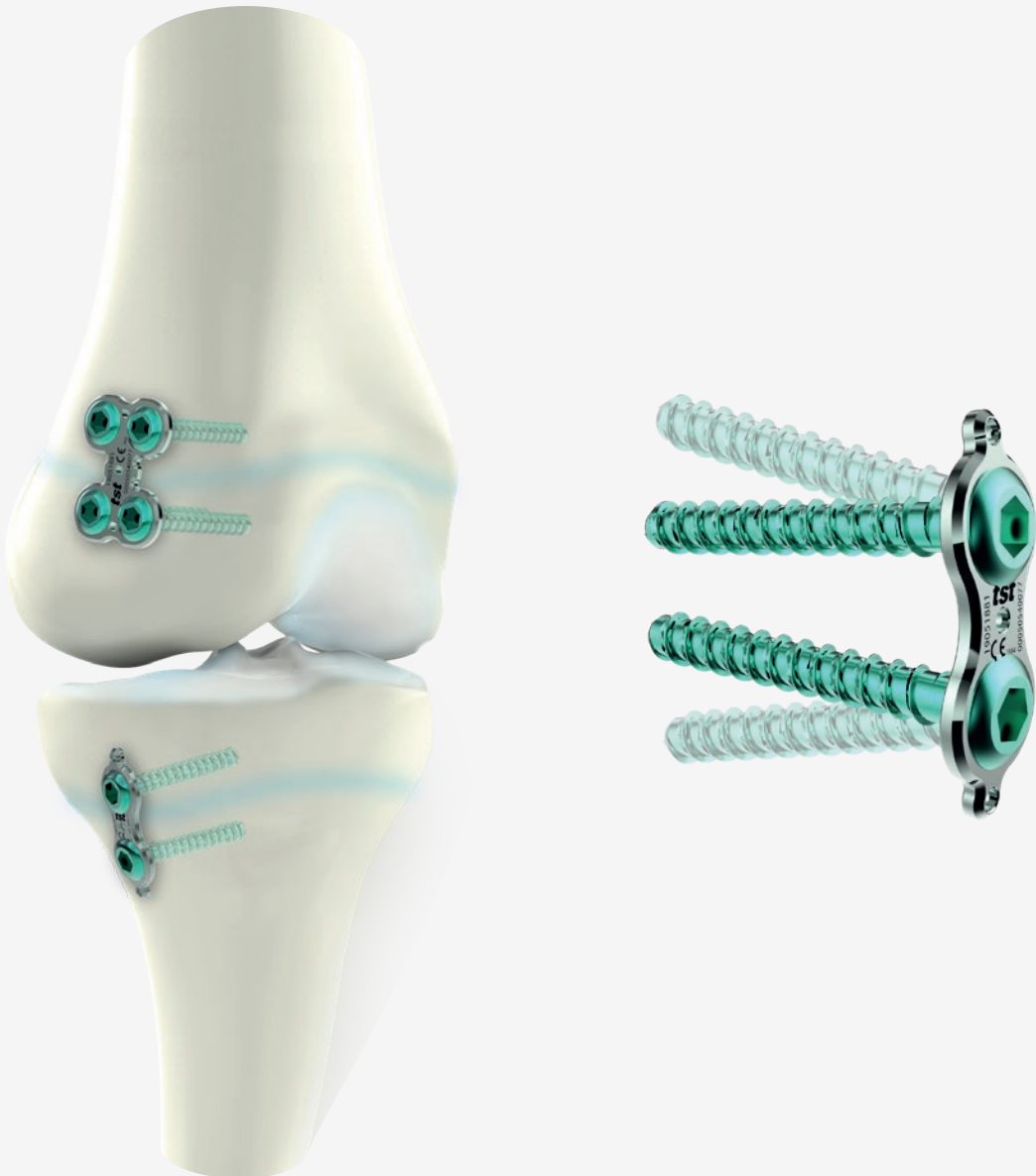

# EIGHT GROWTH PLATES

| TST Code    | Barcode       | Description                | Material | Size  |
|-------------|---------------|----------------------------|----------|-------|
| 37220000012 | 8698673441460 | EIGHT GROWTH PLATE 12 MM   | Ti       | 12 mm |
| 37220000016 | 8698673441477 | EIGHT GROWTH PLATE 16 MM   | Ti       | 16 mm |
| 37222000001 | 8698673469662 | H GROWTH COMPRESSION PLATE | Ti       | 18 mm |
| 37221000015 | 8698673469679 | H GROWTH PLATE 15 MM       | Ti       | 15 mm |



## LOW PROF.SPONG. WIDE HEAD SCREW TI 4.0 MM

| TST Code    | Barcode       | Length | Qty |
|-------------|---------------|--------|-----|
| 20320104014 | 8680858423767 | 14 mm  | 2   |
| 20320104016 | 8680858423774 | 16 mm  | 2   |
| 20320104018 | 8680858423781 | 18 mm  | 2   |
| 20320104020 | 8680858423798 | 20 mm  | 2   |
| 20320104022 | 8680858423804 | 22 mm  | 2   |
| 20320104024 | 8680858423811 | 24 mm  | 2   |
| 20320104026 | 8680858423828 | 26 mm  | 2   |
| 20320104028 | 8680858423835 | 28 mm  | 2   |
| 20320104030 | 8680858423842 | 30 mm  | 2   |
| 20320104035 | 8680858423859 | 35 mm  | 2   |
| 20320104040 | 8680858423866 | 40 mm  | 2   |




## CANN. SPONG. WIDE HEAD SCREW (FULL THRD) 4.0 MM

| TST Code    | Barcode       | Length | Qty |
|-------------|---------------|--------|-----|
| 20324114014 | 8680858427482 | 14 mm  | 2   |
| 20324114016 | 8680858427499 | 16 mm  | 2   |
| 20324114018 | 8680858427505 | 18 mm  | 2   |
| 20324114020 | 8680858427512 | 20 mm  | 2   |
| 20324114022 | 8680858427529 | 22 mm  | 2   |
| 20324114024 | 8680858427536 | 24 mm  | 2   |
| 20324114026 | 8680858427543 | 26 mm  | 2   |
| 20324114028 | 8680858427550 | 28 mm  | 2   |
| 20324114030 | 8680858427567 | 30 mm  | 2   |
| 20324114035 | 8680858427574 | 35 mm  | 2   |
| 20324114040 | 8680858427581 | 40 mm  | 2   |

# Implants & Instruments

| No | TST Code    | Barcode       | Description                                    | Qty |
|----|-------------|---------------|------------------------------------------------|-----|
| 1  | 00000304100 | 8680858437771 | DESIGN TRAY 250X250 MM                         | 1   |
| 2  | 00000020001 | 8680858413089 | EIGHT GROWTH PLATE IMPLANT TRAY                | 1   |
| 3  | 37220000012 | 8698673441460 | EIGHT GROWTH PLATE TI 12 MM                    | 4   |
| 4  | 37220000016 | 8698673441477 | EIGHT GROWTH PLATE TI 16 MM                    | 4   |
| 5  | 37222000001 | 8698673469662 | H GROWTH COMPRESSION PLATE TI                  | 2   |
| 6  | 37221000015 | 8698673469679 | H GROWTH PLATE TI 15 MM                        | 2   |
| 7  | 03020000115 | 8698673493353 | PENSET SMALL                                   | 1   |
| 8  | 03362010010 | 8680858437733 | DRILL SLEEVE HOLDER                            | 1   |
| 9  | 03722000002 | 8699931021967 | EIGHT GROWTH PLATE DRILL SLEEVE                | 1   |
| 10 | 03362007512 | 8680858437726 | K-WIRE SLEEVE 1.2X100 MM                       | 2   |
| 11 | 02050100425 | 8699931005417 | CANNULATED SCREW DRIVER QUICK TIP 2.5XØ5X130MM | 1   |
| 12 | 02020103025 | 8699931010930 | SCREW DRIVER QUICK TIP 2.5XØ5X130 MM           | 1   |
| 13 | 03356200100 | 8680858442256 | LENGTH MEASURING DEVICE FOR K-WIRE 1.2X150 MM  | 1   |
| 14 | 02050100425 | 8698673441187 | SOFT SCREW DRIVER QUICK SMALL                  | 1   |
| 15 | 04551208190 | 8699931029550 | K-WIRE TUBE Ø12XØ8X190 MM                      | 1   |
|    | 23410150112 | 8698673458840 | KIRSCHNER WIRE TROCAR POINT 1.2X150 MM         | 4   |
|    | 23212530170 | 8680858426355 | CANNULATED DRILL 3.2X1.25X30X170 MM            | 1   |
|    | 00270270080 | 8699931009019 | CONTAINER 270X270X80 MM                        | 1   |

# Implants & Instruments

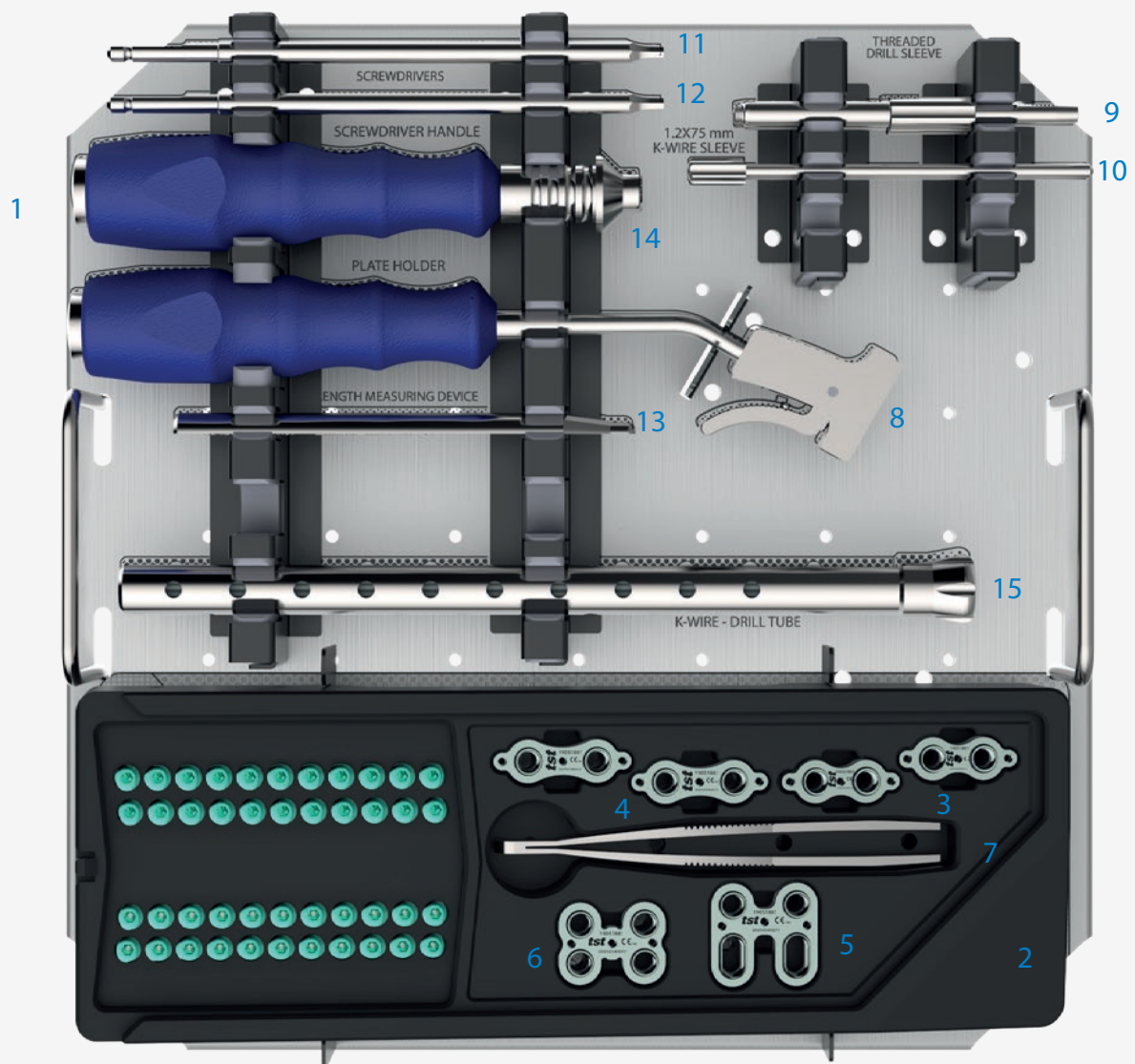

PD-03 Eight Growth Plate Set

