

PD-11
SHELL CALCANEUS
IMPLANT
SET





SHELL CALCANEUS IMPLANT

Planovalgus foot is a major source of disability in children with cerebral palsy (CP). Planovalgus foot is clinically associated with depression of the medial arch, walking disabilities, callosities and ulcerations, and pain. The best one of treatment option in the literature is calcaneal lengthening. Shell Plate was developed for the treatment of planovalgus foot deformity in CP by calcaneal lengthening. The Plate can be inserted into the osteotomy gap based on the amount of talonavicular coverage needed as shown in the picture. If necessary, it can be fixed to the bone with two Locking Peg Screws. Afterward, the osteotomy gap is expanded by turning the inner bolt lateral to the plate.

Indications: Lateral Column Lengthening, Hindfoot Valgus, External Rotation.

TST Code	Barcode	Description	Material	Size
31627010015	8680858441136	SHELL CALCANEUS IMPLANT 15 MM	Ti	15 mm
31627010012	8680858441129	SHELL CALCANEUS IMPLANT 18 MM	Ti	18 mm

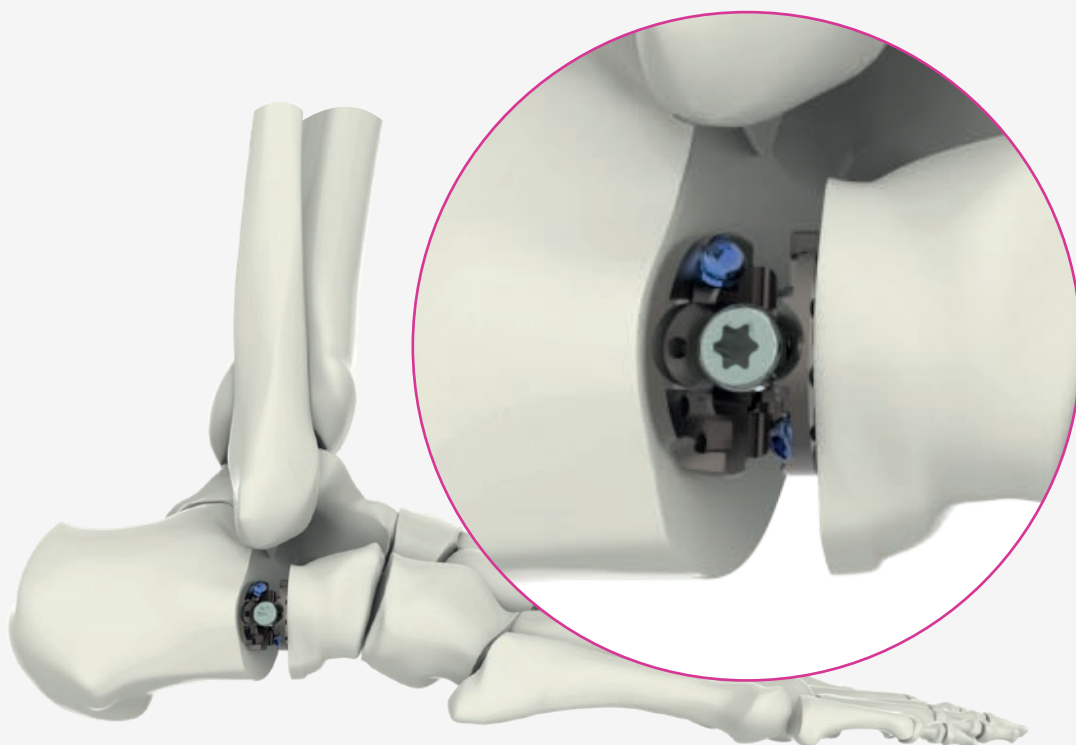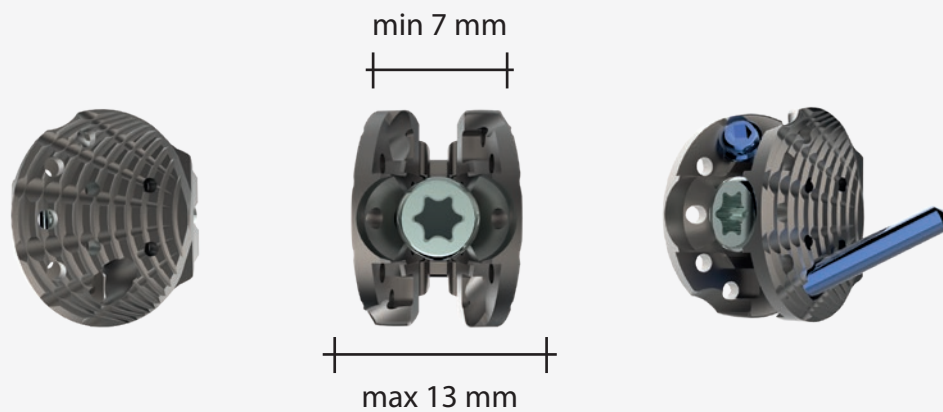

Implants & Instruments

No	TST Code	Barcode	Description	Qty
1	00235235070	8680858451883	CONTAINER 235X235X70 MM	1
2	00000305000	8680858442225	SHELL PLATE & SCREW BOX	1
	31627010015	8680858441136	SHELL CALCANEUS IMPLANT 15 MM	2
	31627010012	8680858441129	SHELL CALCANEUS IMPLANT 18 MM	2
	24627142027	8699931010671	SMALL HEAD LOCK.PEG SCREW TI 2.7X14 MM	4
	24627162027	8699931010688	SMALL HEAD LOCK.PEG SCREW TI 2.7X16 MM	4
	24627182027	8699931010695	SMALL HEAD LOCK.PEG SCREW TI 2.7X18 MM	4
	24627202027	8699931010701	SMALL HEAD LOCK.PEG SCREW TI 2.7X20 MM	4
3	03161700020	8680858441150	SCREWDRIVER STAR TIP T20	1
4	03161700001	8680858441143	SHELL CALCANEUS IMPLANT HOLDER	1
5	03161700002	8680858442249	CALCANEAL SPREADER	1
6	01193003009	8699931028133	BONE HAMMER - SMALL	1
7	02000480017	8699931029253	SCREW DRIVER QUICK HEX. TIP Ø 1.7X80 MM	2
8	02010101001	8698673441187	SOFT SCREW DRIVER QUICK SMALL	1
9	22310130020	8698673450400	BONE DRILL 2.0X130 MM	2
10	02005805020	8680858425730	THREADED DRILL SLEEVE Ø 2X50 MM FOR SMALL HEAD SCR	2

